



Macro163™

*soluble CD163 ELISA assay for the measurement
of macrophage and monocyte activation*

REF²¹ IQP-383

For research use only

▽ 96 tests

📖 package insert

For Research Use Only. Not for Diagnostic Purposes.

©2026 - IQ Products bv. All rights reserved. No parts of these works may be reproduced in any form without permission in writing.

MACRO163™

soluble CD163 ELISA assay for the measurement of macrophage and monocyte activation

Intended use

This Macro163™ kit is intended for the quantification of soluble CD163 (sCD163) in serum or plasma samples. It is designed for research purposes only. The Macro163™ assay has been validated for serum and plasma measurements but can also be applied to other biological fluids like synovial fluid, ascites fluid, cerebrospinal fluid and cell-supernatants. In latter situations adjustments in sample dilutions need to be considered. The Macro163™ assay is completed within 3 hours with less than 1 hour of hands-on-time.

Background information

CD163 is a membrane protein and member of the group B Scavenger Receptor Cysteine-Rich super family specifically expressed on peripheral blood monocytes and macrophages. A particularly high expression is seen in macrophages of the 'alternative activation' phenotype playing a major role in dampening the inflammatory response and in scavenging components of damaged cells. CD163 functions as the receptor for Haptoglobin-hemoglobin complexes, and furthermore CD163 is involved in the regulation of inflammation. Macrophages play a central role in the host response to infection and tissue damage, and are furthermore important in the pathogenesis of autoimmune diseases and cancer. Historically, the presence of CD163 (by immune histochemistry) has been associated with a post-infectious recovery phase and declining inflammation. CD163 positive macrophages constitute a major cell subpopulation in human term placenta suggesting a role for the placenta functioning as an active immunosuppressive biological barrier between mother and fetus.

Importantly, the extracellular part of CD163 is cleaved from the cell-membrane by proteolytic activity upon toll-like receptor activation. This shed CD163 molecule (soluble CD163, sCD163) is detectable in body fluid compartments and is regarded to reflect activation and proliferation of macrophages in inflammatory conditions.

Measurement of sCD163 may be a valuable marker in diseases with macrophage and/or monocyte involvement, such as macrophage activation syndromes (e.g. hemophagocytic syndrome), infections, liver disease, auto-immune disease, atherosclerosis and cancer.

Principle of the Macro163™ assay

The assay is based on the principle of the sandwich ELISA. A polyclonal antibody recognizing CD163 is immobilized on the surface of a microtiterplate. After incubation with the sample or recombinant CD163 standard a second biotinylated monoclonal antibody recognizing CD163 is added. Detection of the latter is done by adding streptavidin-HRP. Using TMB (3, 3', 5, 5'-tetramethylbenzidine) as substrate for the enzyme HRP, the amount of sCD163 protein can be quantified.

Materials supplied with the Macro163™ kit

REF	Item	Description	Amount
A	Reagent A	Recombinant sCD163 standard. Please check the label on the vial for the correct concentration	2x vial
B	Reagent B	Biotinylated detection antibody. Contains 0.02% Sodium azide	135 µL
C	Reagent C	Streptavidine-HRP	150 µL
D	Reagent D (10x)	10x Washing buffer.	40 mL
E	Reagent E (10x)	10x Dilution buffer for samples, standard and controls.	10 mL
F	Reagent F	TMB substrate. Contains <0.02% w/w TMB	11 mL
G	Reagent G	Stop Solution (0.2 M H ₂ SO ₄)	11 mL
H	Pre-coated microtiter plate	12x8 well pre-coated and pre-blocked strips	1x
I	Adhesive plate seals	Plate covers	4x
-	Manual	Instructions for use	1x

Laboratory material and equipment required, but not supplied with the Macro163™ kit

- Laboratory centrifuge for 1.5 mL vials
- 1.5 mL reaction tubes to prepare the recombinant CD163 standard and sample dilutions
- Distilled or demineralized water
- Adjustable micropipettes and tips
- Vortex mixer
- Stopwatch or timer
- Microtiterplate reader with 450 nm (and 620, 650, or 690 nm) filter and software to evaluate ELISA data

Storage

- Reagents stored according to stated storage instructions are stable until the expiration date indicated on the label. Store the reagents for repeated use immediately after usage at 2-8 °C.

Warning and precautions

- For Research Use Only.
- Reagents containing sodium azide (Reagent B [0,002%]) may react with lead or copper plumbing to form explosive metal azides. On disposal, flush with large amounts of water to prevent azide build-up.
- All reagents should be handled in accordance with good laboratory practices using appropriate precautions.
- All samples should be handled with appropriate precautions.
- Do not pipette by mouth and wear gloves during the procedure.
- The test must be performed by well-trained and authorized laboratory technicians.
- Please contact the manufacturer if the original test kit is damaged upon receipt.

Storage of the blood samples

- Blood should be collected in appropriate tubes. Only a few microliters are used in the Macro163™ assay.
- Serum, EDTA-, citrate- or heparin-anticoagulated plasma can be used for analyses.
- There is no interference of EDTA up to 10 mM or concentrations up to 200 mM of hemoglobin.
- Samples should be centrifuged (10 minutes at 1000 g) and serum/plasma drawn off by pipette within 24 hours.
- Soluble CD163 is stable in serum/plasma for more than 2 weeks at 4 °C, over 6 months at -20 °C, and for more than a year at -80 °C.
- After thawing samples can be re-frozen without loss of sCD163 reactivity (7).

Reagent preparation

- Before use, all buffers should be allowed to reach room-temperature (20-25 °C) except for reagent F (TMB) which is kept at 4 °C until section 6.
- Dilute the Washing buffer (reagent D) 10 times with demineralized water before use.
- Dilute the Dilution buffer (reagent E) 10 times with demineralized water before use.

For every Macro163™ 8-well pre-coated strip used we recommend to dilute:

- 0.5 ml reagent E (10x dilution buffer) with demineralized water to a volume of 5 ml,
- 3.5 ml reagent D (10x Washing buffer) with demineralized water to a volume of 35 ml.

Use of a complete Macro163™ pre-coated plate (12 strips) in one assay, we recommend to dilute:

- 6 ml reagent E (10x Dilution buffer) with demineralized water to a volume of 60 ml,
- 40 ml reagent D (10x Washing buffer) with demineralized water to a volume of 400 ml

- **Prior to the Macro163™ assay performance it is recommended to dilute the samples. As the samples may derive from various sources, the following dilutions in 1x Dilution buffer are recommended:**

• Serum/plasma	1:500
• Synovial fluid	1:2500
• Ascites, pleural or pericardial fluid	1:500
• Cerebrospinal fluid	1:50

Macro163™ Assay Procedure

The procedure roughly consists of the following steps

- | | |
|----------------------------------|---------|
| • CD163 standard preparation | step 1 |
| • Sample preparation | step 2 |
| • Sample and Standard incubation | step 3 |
| • Detection-antibody incubation | step 5 |
| • Streptavidine-HRP incubation | step 7 |
| • TMB substrate reaction | step 9 |
| • Adding Stop Solution | step 10 |
| • Data analysis | |

Preparations prior to Macro163™ Assay Procedure

- The incubation steps as described in this section may be performed on a microplate shaker under gentle shaking. This is not necessary for a successful assay performance.
- Prior to their use *reagents A, B and C should be spinned down* for 30 seconds at full speed before opening.

Assay Procedure


1- Preparation of the recombinant sCD163 standard

- Shortly spin the reagent A vial (recombinant sCD163 standard) to make sure all material is at the bottom of the vial before opening. One vial is enough to prepare one standard curve in duplo.
- Prepare reagent A by dissolving the material in 600 µl 1x reagent E (= dilution buffer; referred to as 'DB' in the illustration below). Close the vial and mix thoroughly after addition of the buffer and let rest for 15 minutes.

- **To use the correct concentration of the sCD163 standard in the preparation of the standard dilution curve, check the label on the vial. The concentration after addition of 1x reagent E is stated on the label of the vial.**
- The starting concentration or first point in the standard curve, is the one that is mentioned on the vial (in the illustration noted as 'X').
- From here, the standard dilution curve is a repeated sequence a 2 fold serial dilution described as B – F in the illustration below. 300 µL of the preceding dilution is added to 300 µL of DB. Vortex after mixing.

Note

- For the measurement of samples with very low concentrations of sCD163 (e.g. cerebrospinal fluids or cell culture media) it is advised to discard the value of the highest point standard curve. Instead; make the same dilution curve, but dilute the standard one more step. Start the curve at the second point .

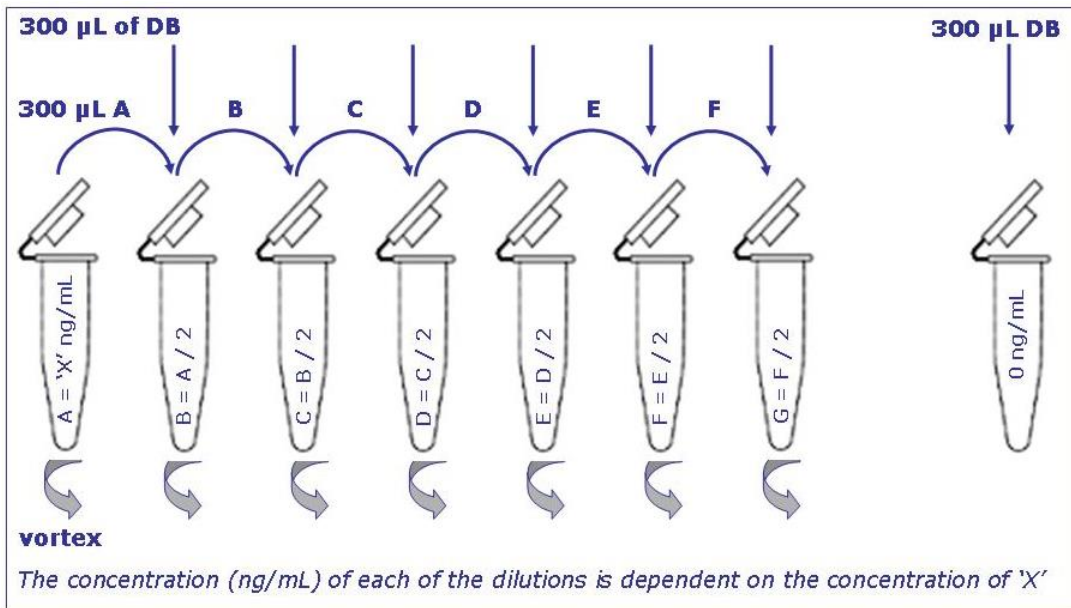

Sample preparation

2– Sample Preparation

- Dilution can be performed by using 1.5 mL reaction tubes.
- Duplicate measurements are recommended.
- When data from various plates will be combined, we recommend to include a set of three to five different samples on each plate, especially when using different lots of kits within one study. By comparing the results of these samples, a correction factor can be determined to make results from different lots comparable.
- Prepare a 1:50 dilution of the serum or plasma sample to be tested by adding 10 µL of sample to 490 µL of Reagent E Dilution buffer. Mix gently by vortexing.
- Prepare a 1:500 final dilution of the sample by adding 50 µL of the 1:50 dilution prepared in step to 450 µL of Reagent E Dilution buffer. Mix gently by vortexing. The volume of the 1:500 sample is sufficient for 4 tests for each sample.

- (This dilution strategy is recommended for serum and plasma samples, resulting in a 1:500 final dilution. Final dilutions can be adjusted to achieve the recommended final dilution for other sample types as described in the table in the Reagent Preparation section.)

3- Adding standards and samples to the pre-coated wells and incubation of samples and standards

- Add 100 μ L of each standard dilution (Reagent A, prepared previously), or sample dilution (prepared in step 2) in duplicate to each well.
- Add 100 μ L of reagent E to two blank wells.
- Mix the fluids gently by tapping the side of the plate or by using a microplate shaker under gentle shaking.
- Seal the plate.
- Incubate at room temperature (20-25C) for 1 hour.
- During this incubation, the detection antibody (reagent B) is prepared (next step).

Preparation of the Detection-antibody (reagent B)

- Dilute the 100x reagent B (100 x detection antibody) in 1x reagent E. For example: Add 100 μ L of 100x reagent B to 10 ml of 1x reagent E. This is sufficient for a full microtiter plate.

4- Washing (manual)

- After the 1 hour incubation the samples/standards are aspirated.
- Invert the plate on a lint-free tissue and tap dry.
- Dispense 200 - 300 μ L 1x reagent D into each well, than aspirate the washing solution.
- Invert plate on a tissue and tap dry.
- Repeat this step for 2 more times.

5- Detection-antibody incubation

- Add 100 μ L of the 1x detection antibody (reagent B) to each well using a multi-channel pipette.
- Seal the plate with the plate sealer.
- Incubate at room temperature (20-25 $^{\circ}$ C) for 1 hour.
- During this incubation the streptavidin-HRP solution (reagent C) is prepared (next step).

Preparation of the streptavidine-HRP solution (reagent C)

- Dilute 100x reagent C in 1x reagent E. For example: Add 100 μ L of reagent C to 10 mL 1x reagent E. This is sufficient for a full microtiter plate.

6- Washing (manual)

- After the 1 hour incubation the diluted reagent B solution is aspirated.
- Invert the plate on a lint-free tissue and tap dry.
- Dispense 200 - 300 μ L 1x reagent D into each well, than aspirate the 1x reagent D.
- Invert plate on a tissue and tap dry.
- Repeat this step for 2 more times.

7- Streptavidine-HRP incubation

- Add 100 μ L 1x streptavidine-HRP solution (reagent C) to each well using a multi-channel pipette.
- Seal the plate with the plate sealer.
- Incubate at room temperature (20-25 $^{\circ}$ C) for 1 hour.
- Aliquot the needed volume of reagent F (number of wells (100 μ l/well) plus an extra 20% of volume) from the stock bottle reagent F with a clean pipette and close the bottle firmly. Wrap it in

aluminum foil and place it in the dark. Reagent F is light sensitive! Place the stock bottle of reagent F back at 2-8 °C.

- Allow the aliquot of TMB to reach room temperature (20-25 °C).

8- Washing (manual)

- After the 1 hour incubation the Streptavidine-HRP solution is aspirated.
- Invert the plate on a lint-free tissue and tap dry.
- Dispense 200 - 300 µL 1x reagent D into each well, than aspirate the 1x reagent D.
- Invert plate on a tissue and tap dry.
- Repeat this step for 4 more times.

9- Adding the TMB


- Add 100 µL of TMB (reagent F) to each well using a multi-channel pipette.
- Seal the plate with a *new* plate sealer.
- Incubate at room temperature (20-25 °C) for 15 min in the dark.

10- Adding the Stop solution

- Add 100 µL Reagent G (Stop Solution) to each well using a multi-channel pipette.
- Mix the fluids gently by tapping the side of the plate or by using a microplate shaker under gentle shaking.
- Read the plate with a 450 nm filter within 15 minutes of adding 0.2 M H₂SO₄ (a 620, 650, or 690 nm filter can be used as reference wavelength).

Data analysis


- Standard curve
The standard curve can be generated by linear interpolation or by automated curve fitting methods as supplied by the standard ELISA software programs. A typical example of a sCD163 ELISA standard curve is shown to the right
- Traceability/correctness
The assay correlates well ($r^2 = 0.97$) with published assays used for clinical studies. The assay was calibrated against a purified human CD163 standard.
- Precision of the Macro163™ assay
Intra-assay variability is 3-6 %; Inter-assay variability is 5-8 % when measuring in duplicates.
- Linearity of the Macro163™ assay
It was shown to be linear from 0.9 – 30 ng/mL.
- Detection level of the Macro163™ assay
Less than 0.23 ng/mL.

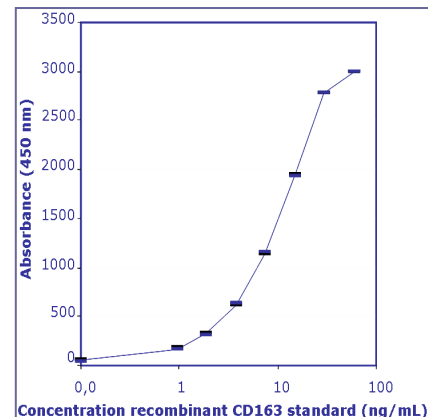

References

- 1) Droste A, et al. Shedding of CD163, a novel regulatory mechanism for a member of the scavenger receptor cysteine-rich family. *Biochem Biophys Res Commun.* 1999 5;256(1):110-3.
- 2) Kristiansen M, et al Identification of the haemoglobin scavenger receptor. *Nature.* 2001 Jan 11;409(6817):198-201.
- 3) Sulahian TH., et al. Development of an ELISA to measure soluble CD163 in biological fluids. *J Immunol Methods.* 2001 Jun 1;252(1-2):25-31
- 4) Møller HJ., et al. Identification of the hemoglobin scavenger receptor/CD163 as a natural soluble protein in plasma. *Blood.* 2002 Jan 1;99(1):378-80.
- 5) Hintz K.A., et al. Endotoxin induces rapid metalloproteinase-mediated shedding followed by up-regulation of the monocyte hemoglobin scavenger receptor CD163. *J Leukoc. Biol.* 2002 Oct;72(4):711-7.
- 6) Møller HJ., et al. Characterization of an enzyme-linked immunosorbent assay for soluble CD163. *Scand J Clin Lab Invest.* 2002;62(4):293-9.
- 7) Møller HJ., et al. Biological variation of soluble CD163. *Scand J Clin Lab Invest.* 2003;63(1):15-21.
- 8) Baeten D, et al. Association of CD163+ macrophages and local production of soluble CD163 with decreased lymphocyte activation in spondylarthropathy synovitis. *Arthritis Rheum.* 2004 May;50(5):1611-23.
- 9) Møller HJ, et al. Plasma level of the macrophage-derived soluble CD163 is increased and positively correlates with severity in Gaucher's disease. *Eur J Haematol.* 2004 Feb;72(2):135-9.
- 10) Moestrup SK, Møller HJ. CD163: a regulated hemoglobin scavenger receptor with a role in the anti-inflammatory response. *Ann Med.* 2004;36(5):347-54. Review. :
- 11) Emmenegger U, et al Haemophagocytic syndromes in adults: current concepts and challenges ahead. *Swiss Med Wkly.* 2005 May 28;135(21-22):299-314. Review.
- 12) Davis BH. And Zarev PV. Human monocyte CD163 expression inversely correlates with soluble CD163 plasma levels. *Cytometry B Clin Cytom.* 2005 Jan;63(1):16-22.
- 13) Fabriek BO., et al. The macrophage scavenger receptor CD163. *Immunobiology.* 2005;210(2-4):153-60. Review.
- 14) Schaer DJ., et al. Soluble hemoglobin-haptoglobin scavenger receptor CD163 as a lineage-specific marker in the reactive hemophagocytic syndrome. *Eur J Haematol.* 2005 Jan;74(1):6-10.
- 15) Aristoteli LP, et al The monocytic lineage specific soluble CD163 is a plasma marker of coronary atherosclerosis. *Atherosclerosis.* 2006 Feb;184(2):342-7. Epub 2005 Jun 23.
- 16) Møller HJ, et al Macrophage serum markers in pneumococcal bacteremia: Prediction of survival by soluble CD163. *Crit Care Med.* 2006 Oct;34(10):2561-6.
- 17) Weiss M, Schneider EM. Soluble CD163: An age-dependent, anti-inflammatory biomarker predicting outcome in sepsis. *Crit Care Med.* 2006 Oct;34(10):2682-3.
- 18) Daly A., et al. Serum levels of soluble CD163 correlate with the inflammatory process in coeliac disease. *Aliment Pharmacol Ther.* 2006 Aug 1;24(3):553-9.
- 19) Bleesing J., et al. The diagnostic significance of soluble CD163 and soluble interleukin-2 receptor alpha-chain in macrophage activation syndrome and untreated new-onset systemic juvenile idiopathic arthritis. *Arthritis Rheum.* 2007 Mar;56(3):965-71.
- 20) Møller HJ, et al. Soluble CD163 from activated macrophages predicts mortality in acute liver failure. *J Hepatol.* 2007 Jun 21; [Epub ahead of print]
- 21) EN ISO 15223-1 Medical devices – Symbols to be used with information to be supplied by the manufacturer - Part 1: General requirements.



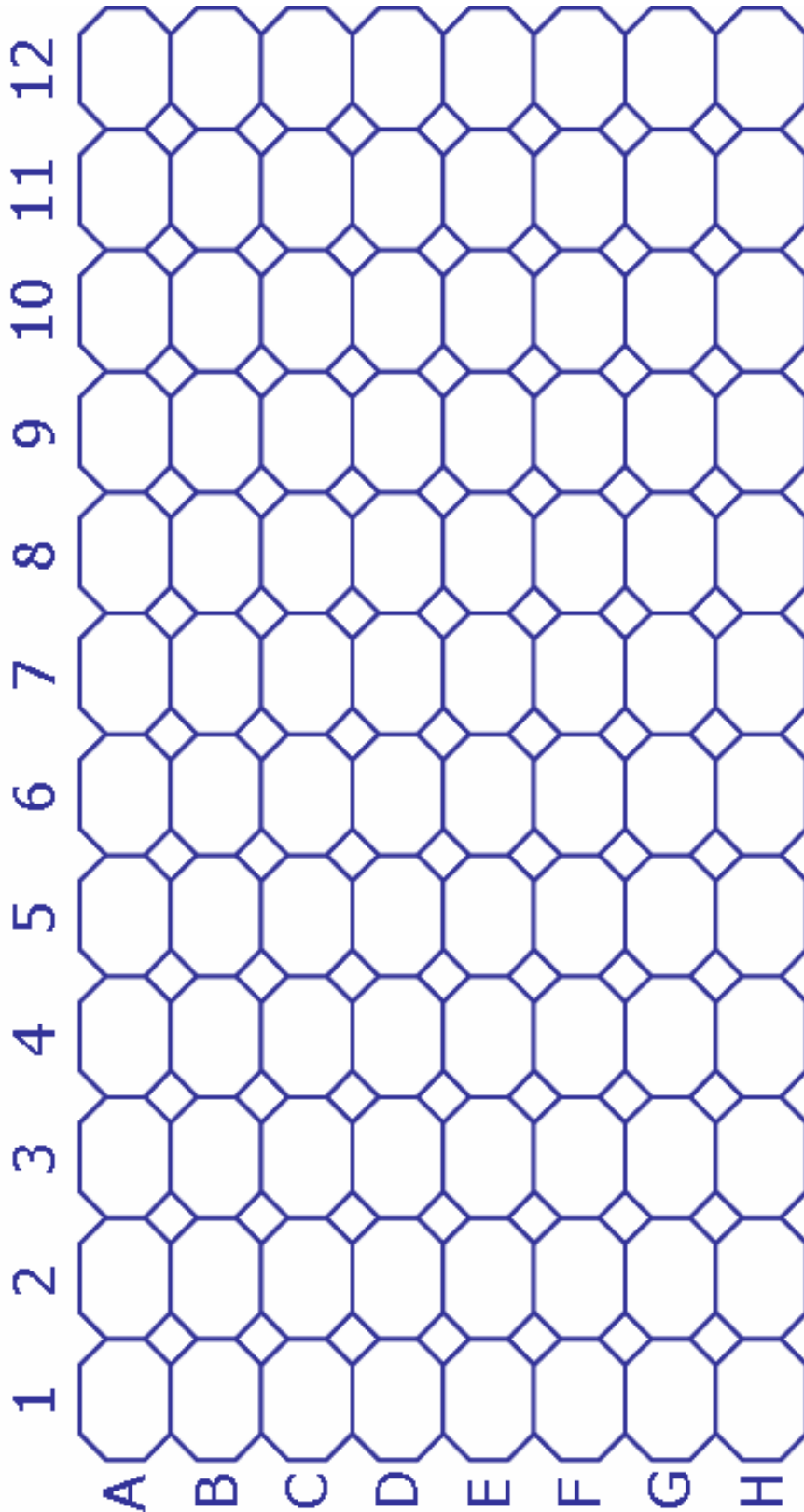







Sampling scheme



Warranty

Products sold hereunder are warranted only to conform to the quantity and contents stated on the label at the time of delivery to the customer. There are no warranties, expressed or implied, which extend beyond the description on the label of the product. IQ Products bv is not liable for property damage, personal injury, or economic loss caused by the product.

Explanation of used symbols

	Consult instructions for use
	Catalogue number
	Sufficient for
	Caution, consult accompanying document
	Keep away from (sun)light
	Biological risks
	Temperature limitation (°C)
	Batch code
	Use by yyyy-mm-dd
	Manufacturer

Contact information

IQ Products BV

www.iqproducts.nl
 Rozenburglaan 13a
 9727 DL Groningen
 The Netherlands
 T 31 (0)50 5757000
 marketing@iqproducts.nl

©2026 - IQ Products bv. All rights reserved. No parts of these works may be reproduced in any form without permissions in writing.