



DOAC-Stop[®]™


The original DOAC eliminator, now also available in micro-tablets and liquid format. Neutralize DOACs from plasma for simplified testing.

DOAC-Stop[®]™ efficiently eliminates DOACs (Direct Oral Anticoagulants)—including dabigatran, apixaban, rivaroxaban and edoxaban—from test plasmas, with minimal effect on currently known clotting variables such as heparin, warfarin, lupus anticoagulants, and clotting factors.



FEATURES & BENEFITS

- Adsorbent product that extracts DOACs from subject plasmas
- Effectively eliminates DOACs from subject plasma for accurate results
- Perform clotting assays on plasma containing DOACs in anticoagulant research
- Identify subject plasmas containing DOACs using simple clotting tests
- Allows for reliable testing of blood coagulation parameters in thrombosis research
- Easy 3-step process: add, mix, centrifuge
- DOAC-Stop[®]™ micro-tablets can be used for samples as small as 0.5ml
- Cost-effective
- DOAC-Stop[®]™ Liquid instantly disperses in samples, allows for customizable volume to meet your sample needs, and does not require centrifugation with mechanical clot-detection instruments
- Current research suggests it can be used with whole blood

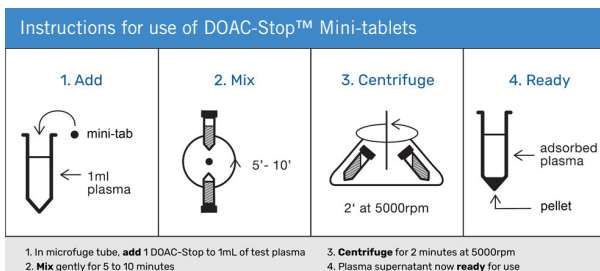

Product	Product #	Description
DOAC-Stop [®] ™ 50 mini-tablets	X9904-50	Each mini-tablet contains enough DOAC-Stop [®] ™ agent for processing 1ml of citrated test plasma (0.5-1.5ml is acceptable). Products are stable at ambient temperature for 3 years if kept dry.
DOAC-Stop [®] ™ 100 mini-tablets	X9904-100	
DOAC-Stop [®] ™ 50 micro-tablets	X9904μ-50	Each micro-tablet contains enough DOAC-Stop [®] ™ agent for processing 0.5ml of citrated test plasma. Products are stable at ambient temperature for 3 years if kept dry.
DOAC-Stop [®] ™ 100 micro-tablets	X9904μ-100	
DOAC-Stop [®] ™ Liquid	X9905-100	DOAC-Stop [®] ™ L is a liquid formulation of the DOAC-Stop [®] ™ mini and micro tabs, made for instant dispersion in samples.
InhibiRef [®]	X9222-D/A/R/E/H/SET	InhibiRef tubes can be used as indicators in clotting and chromogenic assays. Dabi-like, edoxa-like, riva-like, apixa-like and heparin compounds are available individually or as a set of all five.
GO-DOAC Test [®] ™	X9211-FD	GO-DOAC [®] ™ Rapid DOAC Test is a simplified clotting test with high sensitivity for all DOACs, including dabigatran, apixaban, rivaroxaban, and edoxaban.

It's not just support. It's solutions.

800.526.5224 | info@diapharma.com | diapharma.com



*For research use only in US and Canada. Not for use in diagnostic procedures.

ML-00-00349Rev05